

Use & care manual

CORDLESS HANDHELD MASSAGER WITH HEAT

SL-482



 **SNAILAX[®]**



READ ALL INSTRUCTIONS BEFORE USE. SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE.

IMPORTANT SAFEGUARDS:

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

- ALWAYS unplug the appliance from the electrical outlet immediately after using and before cleaning. To disconnect, turn all controls to the 'OFF' position, then remove plug from outlet.
- NEVER leave an appliance unattended when plugged in. Unplug from the mains outlet and turn off the massager when not in use and before putting on or taking off parts or attachments.
- DO NOT reach for an appliance that has fallen in water or other liquids. Switch off at the mains and unplug it immediately. Keep dry – DO NOT operate in wet or moist conditions.
- NEVER insert pins or other metallic fasteners into the appliance.
- Use this appliance for the intended use as described in this booklet. DO NOT use attachments not recommended by SNAILAX.
- NEVER operate the appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to the SNAILAX Service Centre for examination and repair.
- Keep cord away from heated surfaces.
- Excessive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
- NEVER drop or insert any object into any opening.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.
- DO NOT carry this appliance by the cord or use cord as a handle.
- DO NOT attempt to repair the appliance. All servicing of this appliance must be performed by authorised Snailax service personnel only.
- Never block the air openings of the appliance or place it on a soft surface, such as a bed or couch. Keep air openings free of lint, hair etc.
- Do not place or store appliance where it can fall or be pulled into a bath or sink. Do not place in or drop into water or other liquid.
- Failure to follow the above may result in the risk of fire or injury

Features

Cordless Handheld Massager provides vibrating speeds and heat function to meet different massage needs and penetrates deep into your muscle layers to relieve aching and stiff muscles, helping you feel relaxed, recharged and feeling better.

Ergonomic and lightweight design

Control buttons (Intensity, Auto Modes, Heat, Power)

The fixed Round massage head with optional heat function



Note: When you want to activate the heat function, please do not install any of the following interchangeable massage heads

Interchangeable heads for all massage needs

Triangular massage head



Acupuncture massage head



Deep shiatsu massage head



U-shape massage head



Fig.1

TECHNICAL DATA

Model:	SL-482
Voltage:	Input DC 9.0V
Nominal Power:	max. 12 Watt
Battery Max.runtime:	80 minutes

Instructions for use

1. Turn off the massager before replacing another massage head.
2. Your massager comes equipped with 4 attachments (Figure 1). Attach the different heads by twisting clockwise (Figure 2), detach the installed massage head by twisting anti-clockwise and pulling out (Figure 3). Always turn off the device while replacing the heads.

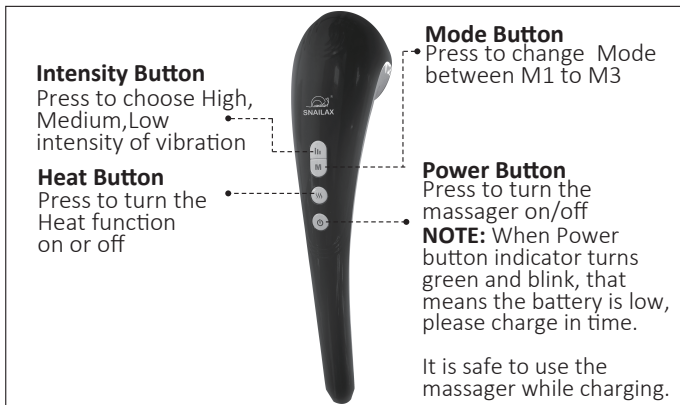


Fig.2



Fig.3

3. Turn on the massager by operating the control buttons (instructions as below)
4. Apply the massager head lightly to the area you wish to treat. Move the unit slowly, holding it over the targeted area. Treat the area for a short period of time and move to another area. NEVER treat the same area for longer than 3 minutes.
5. To turn the unit off by pressing the power button. For your safety the device will automatically turn off after 15 minutes of continuous use.



Battery Charging

1. Connect the charging adapter cable to the corresponding port on the massager



3. Plug the charging adapter into an outlet.
4. Normally it takes 2.5-3 hours to fully charge the battery.
5. Indicator light on POWER Button :
 - When POWER button indicator turns green and blink, that means the battery is low, please charge in time.
 - While charging, the POWER button indicator light will turn 'RED' and blink.
 - When fully charged, the POWER button indicator will keep 'RED' and not blink
 - If the battery is full, the massager can work continuously for around 80 minutes.
6. It is safe to use the massager while charging.

NOTE:

Attention ! Please only use provided charging adapter.

Safety Instructions



READ ALL INSTRUCTIONS BEFORE USE. SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE.

IMPORTANT SAFEGUARDS:

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

- ALWAYS unplug the appliance from the electrical outlet immediately after using and before cleaning. To disconnect, turn all controls to the 'OFF' position, then remove plug from outlet.
- NEVER leave an appliance unattended when plugged in. Unplug from the mains outlet and turn off the massager when not in use and before putting on or taking off parts or attachments.
- DO NOT reach for an appliance that has fallen in water or other liquids. Switch off at the mains and unplug it immediately. Keep dry – DO NOT operate in wet or moist conditions.
- NEVER insert pins or other metallic fasteners into the appliance.
- Use this appliance for the intended use as described in this booklet. DO NOT use attachments not recommended by SNAILAX.
- NEVER operate the appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to the SNAILAX Service Centre for examination and repair.
- Keep cord away from heated surfaces.
- Excessive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
- NEVER drop or insert any object into any opening.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.
- DO NOT carry this appliance by the cord or use cord as a handle.
- DO NOT attempt to repair the appliance. All servicing of this appliance must be performed by authorised Snailax service personnel only.
- Never block the air openings of the appliance or place it on a soft surface, such as a bed or couch. Keep air openings free of lint, hair etc.
- Do not place or store appliance where it can fall or be pulled into a bath or sink. Do not place in or drop into water or other liquid.
- Failure to follow the above may result in the risk of fire or injury.
- Do not use this product for medical treatment.

Safety Cautions

PLEASE READ THIS SECTION CAREFULLY BEFORE USING THE APPLIANCE.

If you have any concerns regarding your health, consult a doctor before using this appliance.

- Individuals with pacemakers and pregnant women should consult a doctor before using this appliance. Not recommended for use by Diabetics
- DO NOT use on an infant, invalid or on a sleeping or unconscious person.
- DO NOT use on insensitive skin or on a person with poor blood circulation.
- This appliance should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls or who has sensory deficiencies in the lower half of their body.
- NEVER use directly on swollen or inflamed areas or skin eruptions.
- If you feel any discomfort whilst using this appliance, discontinue use and consult your doctor.
- This is a non-professional appliance, designed for personal use and intended to provide a soothing massage to worn muscles. DO NOT use as a substitute for medical attention.
- Do not use for longer than the recommended time.
- If the product feels excessively hot, switch off at mains immediately.
- You may soften massage force by placing a towel between yourself and the unit.
- Please ensure that all hair, clothing and jewellery are kept clear of the massage mechanism or any other moving parts of the product at all times.

Care and Maintenance

- For maintenance and repair of this product, please contact our service center.
- Please turn off the device immediately when not in use.
- Place this product in a dry, dust-free place.
- Do not store this product near high temperature or open flame.
- Do not expose it to direct sunlight, it will cause discoloration and damage.
- Avoid the use of benzene, thinners, or any other corrosive solvents, or liquid on this product and its accessories.
- Always unplug from electrical outlet and turn off the device before cleaning this product and its accessories,
- Wipe the product with a soft, dry cloth to clean the machine.

TROUBLESHOOTING

FAULT	SOLUTION
THE DEVICE FAIL TO START	CHECK IF THE BATTERY IS OUT OF POWER
	CHECK IF THE MODE OR HEAT BUTTON IS TURNED ON.
THE VIBRATION BECOMES WEAK	THE BATTERY IS LOW, PLEASE CHARGE IN TIME
THE DEVICE SWITCHES OFF SUDDENLY	BATTERY POWER RUN OUT, PLEASE RECHARGE IN TIME

Warranty

If you have any issue about the product, please feel free to contact us by sending email to support@snailax.com We will strive to provide the best service possible within 24 hours.

30 days unconditionally Return

Snailax product can be returned to receive a full refund for any reason within 30 days. Please contact our customers service (support@snailax.com), our staff will contact you within 24 hours.

90 days refund/replace

Snailax product can be returned / replaced within 90 days if the product break down in the period of proper use.

12 months Warranty

If the product break down within 12 months in the period of proper use, customers can enjoy the relevant warranty to get them replaced.

Attention!

No warranty will be given to any force majeure and man-made causes to a defective product, such as improper care, personal tear down and willful damage, etc.

Extend Warranty for Free

1) Enter the following URL or scan the QR code below to find Snailax facebook page and like it, enter“Warranty” to messenger to extend your warranty from 1 year to 3 years.



<https://www.facebook.com/snailax>

OR 2) Send message“Warranty” and email us support@snailax.com to extend your warranty from 1 year to 3 years.

Have question ?



Tel: 734-709-6982

Monday-Friday 9:00AM-4:30PM



Email: support@snailax.com

FCC STATEMENT

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio / TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this unit is not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Add:2750 Carpenter Road, Suite 1 B, Ann Arbor, 48108

Tel:734-709-6982

support@snailax.com

www.snailax.com