

Use & care manual

HAND MASSAGER

SL-488



 **SNAILAX[®]**

Safety Instructions



READ ALL INSTRUCTIONS BEFORE USE. SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE.

- ALWAYS unplug the appliance from the electrical outlet immediately after using and before cleaning.
- Do not wet, do not use pins, never remove cover.
- This item is NOT A TOY. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.
- This appliance should never be left unattended when plugged in.
- Any possible repairs may only be carried out by authorized expert staff. Improper use and unauthorized repairs are not permitted for safety reasons and lead to loss of warranty.
- Never touch the power plug with wet hands.
- Please avoid contact of the device with water, high temperatures and direct sunlight.
- Do not use any damaged cables, plugs or loose sockets.
- In case of malfunction, disconnect immediately from the mains.
- Do not use if you have skin disorders, open wounds, or swollen or inflamed areas.
- Misuse or incorrect use excludes any liability for damage.
- Do not use this device while driving.
- Do not use while you sleep.
- Do not use this product for medical treatment.
- In order to avoid excessively stimulating the muscles and nerves, the recommended massage time should not be exceeded 15 minutes at a time.

If you suffer from any of the medical conditions listed below, consult with a doctor before using this relaxation massage device.

If you:

suffer from a heart condition

have a sensory disorder or pacemaker

suffer from osteoporosis or diabetes

have a malignant tumor

have neurological disorder or perceptual disorder on head

have fractured or head surgery of your massaging area

head skin disease and allergic dermatitis

suffering from acute diseases or are undergoing medical treatment

we advise you to consult with your doctor before using this massager.

- The packaging material supplied should not be used as a toy.

Failure to observe the above instructions may constitute misuse of the product and can cause serious injury or burns

Features

Thank you for purchasing SNAILAX Hand Massager

With normal care and proper treatment, it will provide years of reliable service

FEATURES

1. Compression Massage with soothing heat function.
2. Six massage modes to choose
3. Each massage mode with 6 intensity levels for options
4. Build in rechargeable Lithium battery
5. Optional vibration function

CONTENTS

- Hand Massager
- USB Charging cable



TECHNICAL DATA

Model:	SL-488
Weight:	2.13 Lbs
Adapter Input:	100-240V~ 50/60 Hz
Rated Voltage:	DC5V 1A
Rated Power:	5W
Default Automatic Runtime:	15 minutes
Battery Spec:	3.7V Lithium battery, 2500mAh

Intensity button


Press to cycle intensity from L3- L4- L5- L6- L1- L2, and the corresponding indicator will light on. Default intensity is L3.

Heat button

Press to turn ON/OFF heat function

Power button

Press to turn on/off the device.

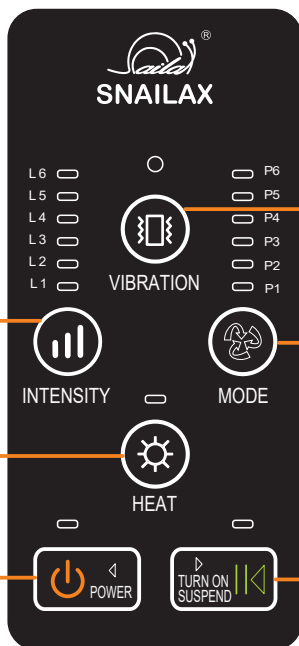

Vibration button

Press to turn on vibration function, you may cycle through low, high or off setting, the indicator light will turn green/blue/off respectively.

Mode button

Press to select message mode between P1 to P6, the corresponding indicator will light on. Default mode is P1

Suspend button

Press to pause the massager, the intensity and mode indicators will begin to flash, press again to return to the previous working state

! Note

This device will switch off automatically after 15minutes timer run out. You should allow your body some rest periods. To avoid straining your muscles too much, we recommend to not exceed a continuous massage of **15** Minutes per time and not exceed **3** times use one day

Battery Charging

1. When the battery is low, all indicators will start flashing.
2. Connect the USB cable to any USB charger, it can be plugged in the portable power bank, mobile phone battery chargers, laptop and any USB port.



3. When charging, the power indicator light will start flashing.
4. When full charge, the power indicator turn off.
5. Normally it takes about 5 hours to fully charge.
6. Machine will be turned off automatically when charging.
7. After full charged, the device can work continuously for around 1.5 hours (heating function is on). or around 4.5 hours (heating function is off)

Note: before using this product for the first time, please fully charge it. If the device is not in use for a long time, the product should be charged once a month to protect the battery life from shortening.

Battery Disposal Warning:

The unit contains a Lithium Ion battery.

Please contact your local recycling center to find out how to properly dispose of the product.

Troubleshooting

Trouble	Reason/Solution
Can not start the massager	Battery Low, Please charge it before using
The massager stoped abnormally	Battery runs out of power, charge the battery; Self-timed 15 mins is out, restart the product
Massager slight noise	It is normal working noise of the internal mechanism

Care and Maintenance

Care and cleaning instructions

1. Always unplug the unit before cleaning it.
2. Before cleaning, allow unit to cool.
3. Dampen a cloth in water or a 3%–5% mild detergent solution.
4. Wipe the dirty areas with the wet cloth.
5. Do not use sharp brush, sharp tool, gasoline and diesel to clean it.
6. Do not use wet cloth clean control panel, lest cause failure, only use dry cloth to wipe the control panel.
7. Never submerge any part of the unit in liquid.
8. Wait for the unit to become fully dry before use.

Maintenance

1. Please keep it into a box and keep in a safe, dry and cool place.
2. Do not keep it in any place with high temperature, near fire or direct sunlight, or it may faded or damaged
3. Do not contact benzene, dilient or any other corrosive liquid.

Warranty

If you have any issue about the product, please feel free to contact us by sending email to support@snailax.com. We will strive to provide the best service possible within 24 hours.

30 days unconditionally Return

Snailax product can be returned to receive a full refund for any reason within 30 days. Please contact our customers service (support@snailax.com), our staff will contact you within 24 hours.

90 days refund/replace

Snailax product can be returned / replaced within 90 days if the product break down in the period of proper use.

12 months Warranty

If the product break down within 12 months in the period of proper use, customers can enjoy the relevant warranty to get them replaced.

Attention!

No warranty will be given to any force majeure and man-made causes to a defective product, such as improper care, personal tear down and willful damage, etc.

Extend Warranty for Free

1) Enter the following URL or scan the QR code below to find Snailax facebook page and like it, enter "Warranty" to messenger to extend your warranty from 1 year to 3 years.



<https://www.facebook.com/snailax>

OR 2) Send message "Warranty" and email us support@snailax.com to extend your warranty from 1 year to 3 years.

Have question ?



Tel: 734-709-6982

Monday-Friday 9:00AM-4:30PM



Email: support@snailax.com

FCC STATEMENT

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
 - Increase the separation between the equipment and receiver.
 - Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
 - Consult the dealer or an experienced radio / TV technician for help.
- This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this unit is not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.



2750 Carpenter Road,
Suite 1 B, Ann Arbor, 48108
Tel: 734-709-6982
support@snailax.com
www.snailax.com