Angelcare®

IMPORTED BY / IMPORTÉ PAR:
ANGELCARE CANADA INC.
2000, MCGILL COLLEGE AVENUE, SUITE 250
MONTREAL, QUEBEC, H3A 3H3, CANADA
WWW.ANGELCAREBABY.COM
MADE IN CHINA / FABRIQUÉ EN CHINE

IMPORTANT! READ CAREFULLY AND KEEP FOR FUTURE REFERENCE

Angelcare Baby Bath & Room Thermometer

The Angelcare Baby Bath & Room Thermometer allows you to easily determine the ideal temperature for your baby's bath water or room. Gently shake the thermometer to turn it ON, the temperature will be displayed and it is ready to use. Without any activity, after 30 minutes, it will automatically turn OFF to power saving mode.

For use as a Bath Thermometer

Fill the bath to a safe level for your child and ensure the water is thoroughly mixed to avoid hot areas and to obtain an accurate reading. Place the thermometer face up in the bath water. The temperature will be displayed within a few seconds.

At a temperature of 39.0°C or higher the LCD screen will display 'HOT' and a red light will flash to alert that the water temperature is too high for a baby and could get burnt.

A comfortable bath temperature for a baby is between 35.0°C and 38.0°C.

At a temperature of 32.0°C or lower the LCD screen will display 'COLD' to alert that the water temperature is too low for a baby.

The thermometer is safe for your baby to play with in the bath with parental supervision as it is designed to Toy Safety Standards.

On removal of thermometer from the bath, clean and wipe with a dry cloth or towel and only use non-abrasive cleaning agents such as soapy water if required. Do not immerse in water when cleaning.

For use as a Room Thermometer

Gently shake to activate.

19°C is the ideal room temperature for a baby when sleeping.

The word 'COLD' appearing at the top of the display screen is only applicable when the thermometer is used to measure the temperature of the bath water.

! WARNING

This product contains a Button or Coin Cell Battery. A swallowed Button or Coin Cell Battery can cause internal chemical burns in as little as two hours and lead to death. Dispose of used/waste batteries immediately. Keep new and used batteries away from children. If you think batteries might have been swallowed or placed inside any part of the body, seek immediate medical attention.

CAUTION - Do not leave the child without adult supervision.

Always check temperature of the water before placing baby in the bath.

Only to be used in water in its fully assembled state.

Remove the battery during long periods of non-use.

Leakage and corrosion can damage the product.

Use a battery of the same or equivalent type and ensure correct polarity

when inserted noting the positive terminal (+) always face up in battery compartment.

Do not short circuit battery terminals.

Always make sure battery compartment is secure and sealed.

Do not use rechargeable batteries.

Changing the Battery

After about a year the CR2032 battery may go flat and the display may dim or not appear.

Open the cover at the bottom of the thermometer using a coin or similar

in the direction indicated on the cover.

Insert the battery and ensure the positive sign on the battery (+)

is facing up (diagram Fig.1).

Replace cover and ensure the seal is in place and tighten cover

in the opposite direction to opening as indicated on the cover.

If the seal is damaged in any way or letting in water stop using the product immediately.

Technical Data

Model BT-01 | 1 x CR2032 DC 3.0V Battery | Temp. accuracy: ±1°C | Temp. measurement and storage range: 0°C - 50°C, no restriction in humidity. | Display resolution: 0.1°C | IP27 rated

PROTECT THE ENVIRONMENT - The wheelie bin symbols mean that according to local laws and regulations your product and/ or its battery shall be disposed of separately from household waste. When the product/battery reaches its end of life take it to a collection point designated by local authorities. The separate disposal and recycling will help conserve natural resources and protect human health and the environment.

Electromagnetic Fields (EMF) - The Angelcare Baby Bath & Room Thermometer complies with all standards regarding electromagnetic fields (EMF). If handled correctly and in accordance to user instructions supplied, the device is safe to use based on scientific evidence available today.

Contact information

IMPORTANT! This Angelcare Baby Bath & Room Thermometer is a personal care product. Please do not return this product to the store. If you have questions or need help, please call your distributor before returning this product.

We, at Angelcare®, stand behind all of our products. If you are not completely satisfied or have any questions, please contact your distributor before returning this product.

United States and Canada: Please contact us at: info@angelcarebaby.com

Other countries: Find your local distributor by visiting our website: angelcarebaby.com/distributors

To register your product, visit our online registration at: www.angelcarebaby.com.

We will use the information provided to contact you if there is a Safety alert or Recall for this product.

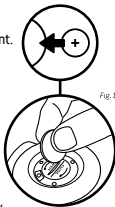
Fig.1

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