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
area of the body that is swollen, burned, bruised, inflamed or where skin eruptions, ulcerated sores or varicose veins are present.

- This product should never be used by any individual suffering from any ailments that would limit or altogether inhibit the user's capacity to feel or have any sensation in the body part that this product is being used for.
- Use of this product should be pleasant and comfortable. Should pain or discomfort result, discontinue use and consult your physician.
- Ensure that your hands are dry when operating the controls or removing the plug.

HOW TO ASSEMBLE

To assemble the heating pad, insert the electrical cord into the connecting port on the side of the heating pad. If moist heat is desired, follow instructions in the Moist Heat Treatment section.

HEATING PAD FEATURES

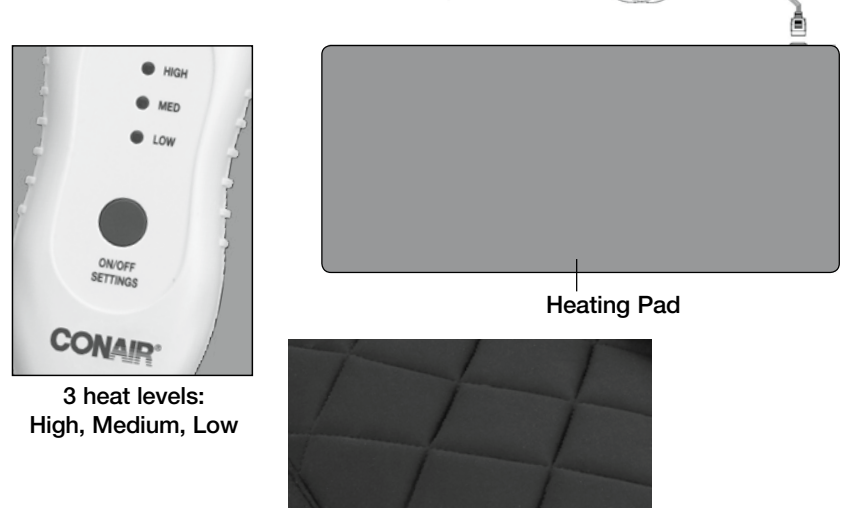MACHINE WASHING INSTRUCTIONS

Wash by hand or as follows in an automatic washer using slow agitation and normal spin speed: Fill washer with warm water. Add minimum amount of all-purpose detergent, agitate to dissolve, measure heating pad before washing. Add heating pad and soak 5 minutes. Agitate 2 minutes. Spin. Fill for cool rinse, agitate 1 minute, then spin.

DRYING INSTRUCTIONS

NEVER use commercial dryers at the laundromat, as they may be extra hot and could damage the wiring of the product. DO NOT iron. DO NOT plug in until completely dry. To machine dry, place in automatic dryer on high heat for only 5 minutes then line dry (recommended). To completely dry, place item in a dryer on low heat and allow to tumble until dry. Ensure that the textile has enough room to tumble freely. If product draws during washing, gently stretch it while it is still damp until it returns to its original size. DO NOT iron your heating pad. To line dry, simply drape the item over two parallel clotheslines or a shower rod. DO NOT reconnect the item until it is completely dry. DO NOT use clothespins, as they can pinch the wires.

STORAGE

Unplug and store this item as you would any other heating pad, folded, and in a dry place.

POLARIZED PLUG

This appliance has a polarized plug (one blade is wider than the other) as a safety feature. To reduce the risk of electric shock, this plug is intended to fit into a polarized outlet only one way. If the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician. Do not attempt to defeat this safety feature. If the plug fits loosely into the AC outlet or if the AC outlet feels warm, do not use that outlet.

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
2. Do not touch the plug or device with wet hands.
3. Do not touch unit to wet skin.
4. Do not use if cord or plug appears to be damaged.
5. Do not use the unit if it has been damaged, torn, or has come in contact with water.
6. Do not pull, lift, or carry unit by the power cord.
7. Do not operate the unit in areas where flammable or combustible products or vapors may be present. Explosion or fire may occur.
8. Do not allow children to use – this is not a toy. Do not leave small children unattended close to unit.
9. Not to be used by disabled, sleeping or unconscious persons.

Cautions:

- This heating pad should not be used in the following cases: If you have diabetes, tuberculosis, benign or malignant tumors, phlebitis and thrombosis, hemorrhages, inflammation of the skin, undetermined calf pain or joint pain, or if you have poor circulation. Consult a physician before using if you have any kind of fungal infection, swelling, fracture or persistent pain.
- Consult your physician should you have any therapeutic questions regarding use of this heating pad.
- Do not use on any unexplained pain, or swollen muscles, or following a serious muscular injury before consulting your physician.
- Keep away from children.
- Discontinue use of product if pain increases or swelling occurs after application.
- Do not use while bathing or showering. Use only in a dry area.
- Consult your doctor before using in case of illness or medical condition, including but not limited to pacemaker use, pregnancy or cancer.
- If you experience pain in a muscle or joint for a prolonged period of time, discontinue use and consult your doctor. Persistent pain could be a symptom of a more serious condition.
- Never use any heating products on open wounds, discolored areas, or any

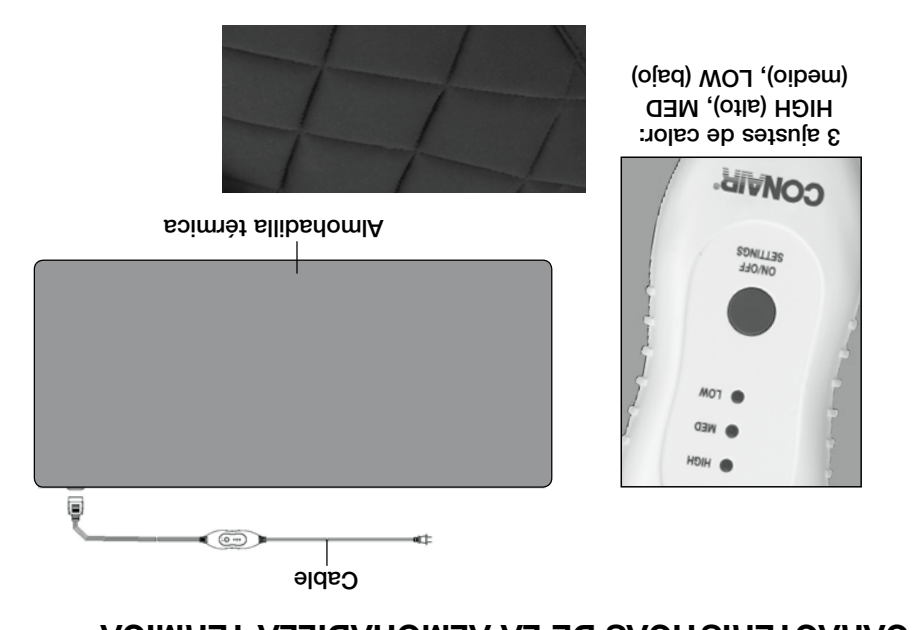INSTRUCCIONES DE LAVADO A MÁQUINA

USE, DISCARD THE PAD IF INNER COVERING OR LINE CORD SHOW ANY SIGN OF DETERIORATION.


15. DO NOT TAMPER WITH THIS PAD IN ANY WAY. THERE ARE NO USER SERVICEABLE PARTS. IF, FOR ANY REASON, THIS PAD DOES NOT FUNCTION SATISFACTORILY, RETURN IT TO AUTHORIZED SERVICE CENTER AT ADDRESS GIVEN.
16. DO NOT USE THIS PAD WITH LINIMENT, SALVE OR OINTMENT PREPARATIONS THAT CONTAIN HEAT-PRODUCING INGREDIENTS, AS SKIN BURNS COULD RESULT. USE WHEREVER HOT APPLICATIONS ARE DESIRABLE FOR PERSONAL COMFORT AND WHEREVER RECOMMENDED BY YOUR PHYSICIAN OR DENTIST FOR THE RELIEF OF PAIN.
17. LOOP CORD LOOSELY WHEN STORING; TIGHT WRAPPING MAY DAMAGE CORD AND INTERNAL PARTS.
18. KEEP THE PAD AWAY FROM HEATED SURFACES.
19. PLACE PAD ON TOP AND NOT UNDER THE PART OF THE BODY NEEDING HEAT.
20. DO NOT PUT THIS PAD IN THE MICROWAVE.
21. DO NOT LEAVE PAD UNATTENDED, ESPECIALLY IF CHILDREN ARE PRESENT.
22. DO NOT USE PADS UNDERNEATH BEDCOVERS AND OTHER TYPES OF BLANKETS. DO NOT PLACE COVERS, SUCH AS SHEETS AND BLANKETS, ABOVE PADS. ONLY THE CLOTH COVER PROVIDED WITH THE PAD SHOULD BE USED.
23. DO NOT USE PAD ON SENSITIVE SKIN OR WHEN TAKING MEDICATION THAT CAN CAUSE SENSITIVE SKIN.
24. DO NOT DRY-CLEAN THIS HEATING PAD. CLEANING SOLVENTS MAY HAVE A DETERIORATING EFFECT ON THE INSULATION OF THE HEATING ELEMENT.

WARNING – The Conair Memory Foam Heating Pad is an electrical device. To avoid electric shock:

1. Do not wash the unit with water. See Washing Instructions.

IMPORTANT SAFETY INSTRUCTIONS

- DANGER** – To reduce the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions:
1. READ ALL INSTRUCTIONS.
 2. ALWAYS UNPLUG THE UNIT FROM THE ELECTRICAL OUTLET IMMEDIATELY AFTER USING.
 3. DO NOT USE WHILE SLEEPING. AVOID FALLING ASLEEP WHILE USING THE PRODUCT.
 4. DO NOT USE ON AN INFANT OR AN ANIMAL.
 5. THIS PAD IS NOT TO BE USED ON OR BY ANY INDIVIDUAL WITH CERTAIN DISABILITIES, SLEEPING OR UNCONSCIOUS PERSON, OR A PERSON WITH POOR BLOOD CIRCULATION UNLESS CAREFULLY ATTENDED.
 6. DO NOT USE ON AREAS OF INSENSITIVE SKIN.
 7. BURNS CAN OCCUR REGARDLESS OF CONTROL SETTING; CHECK SKIN UNDER PAD FREQUENTLY.
 8. NEVER USE PAD WITHOUT REMOVABLE COVER IN PLACE.
 9. DO NOT USE IN AREAS WHERE OXYGEN IS BEING ADMINISTERED.
 10. DO NOT USE PINS OR OTHER METALLIC MEANS TO FASTEN THIS PAD IN PLACE.
 11. DO NOT SIT ON OR CRUSH PAD – AVOID SHARP FOLDS.
 12. NEVER PULL THIS PAD BY THE SUPPLY CORD.
 13. DO NOT USE THE CORD AS A HANDLE.
 14. CAREFULLY EXAMINE INNER COVER AND LINE CORD BEFORE EACH

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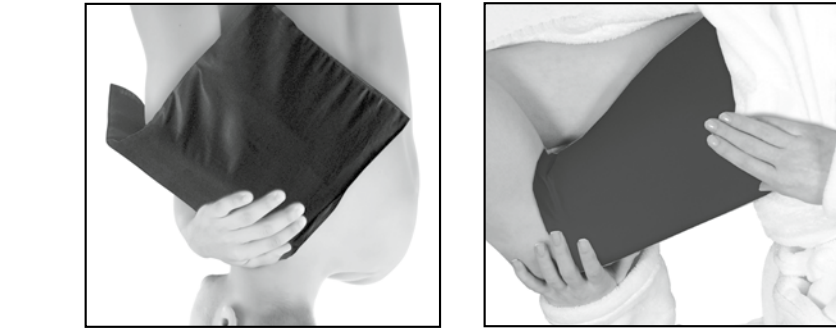Do not use under covers

Do not sit on heating pad

Do not sharply fold or crease heating pad

Not for use by children

Do not sleep while using heating pad

Do not use heating pad punched up

For additional warnings, please see the Important Safety Instructions and Caution sections.

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Centro de servicio
Conair Corporation
Service Department
7475 North Glen Harbor Blvd.
Glendale, AZ 85307

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Conair reparará o reemplazará (a su opción) your unit free of charge for 36 months from the date of purchase if the appliance is defective in workmanship or materials.

To obtain service under this warranty, return the defective product to the service center listed, together with your purchase receipt. California residents need only provide proof of purchase and should call 1-800-3-CONAIR for shipping instructions. In the absence of a purchase receipt, the warranty period shall be 36 months from the date of manufacture.

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©2020 Conair Corporation