# Shiatsu Massage Chair Traveling Shiatsu Mechanism

## Instruction Manual and Warranty Information

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**2** year  
limited warranty

## IMPORTANT SAFEGUARDS

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING:

READ ALL INSTRUCTIONS BEFORE USING

### **DANGER** - TO REDUCE THE RISK OF ELECTRIC SHOCK:

- ALWAYS unplug the appliance from the electrical outlet immediately after using and before cleaning.
- DO NOT reach for an appliance that has fallen into water. Unplug it immediately.
- DO NOT use while bathing or showering.

### **WARNING** - TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK, FIRE OR INJURY TO PERSONS:

- An appliance should never be left unattended when plugged in. Unplug from outlet when not in use and before putting on or taking off parts or attachments.
- Close supervision is necessary when this appliance is used by, on or near children, invalids or disabled persons.
- Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by HoMedics; specifically any attachments not provided with the unit.
- NEVER operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return it to HoMedics Service Center for examination and repair.
- Keep cord away from heated surfaces.
- NEVER use while sleeping or fall asleep while using the massager.
- NEVER drop or insert any object into any opening.
- DO NOT operate where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate under a blanket or pillow. Excessive heating can occur and cause fire, electrocution or injury to persons.
- DO NOT carry this appliance by supply cord or use cord as handle.

- To disconnect, turn all controls to the "OFF" position, then remove plug from outlet.
- Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- Do not use outdoors.

## SAVE THESE INSTRUCTIONS

### **Caution— Please read all instructions carefully before operating.**

- This product is not intended for medical use. It is intended only to provide a luxurious massage.
- Consult your doctor prior to using this product, if
  - You are pregnant
  - You have a pacemaker
  - You have any concerns regarding your health
- Not recommended for use by diabetics.
- NEVER leave the appliance unattended, especially if children are present.
- NEVER cover the appliance when it is in operation.
- DO NOT use this product for more than 15 minutes at a time.
- Extensive use could lead to the products excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
- NEVER use this product directly on swollen or inflamed areas or skin eruptions.
- DO NOT use this product as a substitute for medical attention.
- DO NOT use this product before bed. The massage has a stimulating effect and can delay sleep.
- This product should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls.
- This product should never be used by any individual suffering from any ailments that would limit or altogether inhibit the user's capacity to feel or have any sensation in the body part that the product is being used on.
- This unit should not be used by children or invalids without adult supervision.


## Instructions for Use

To reduce the risk of electric shock, this appliance has a polarized plug (one blade is wider than the other). This plug will fit in a polarized outlet only one way. If the plug does not fit fully in the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install the proper outlet. Do not change the plug in any way.

# Assembly Instructions



## INVENTORY OF PARTS:

- |                           |                             |
|---------------------------|-----------------------------|
| A. Seat Back (x1)         | G. Seat Base Mechanism (x1) |
| B. 6mm Allen Wrench (x1)  | H. Plastic Cone Cover (x1)  |
| C. Armrests (x2)          | I. Star Base (x1)           |
| D. Screw Covers (x8)      | J. Rolling Castors (x5)     |
| E. 8m x 25mm Screws (x12) | K. Pneumatic Cylinder (x1)  |
| F. Seat Cushion (x1)      |                             |

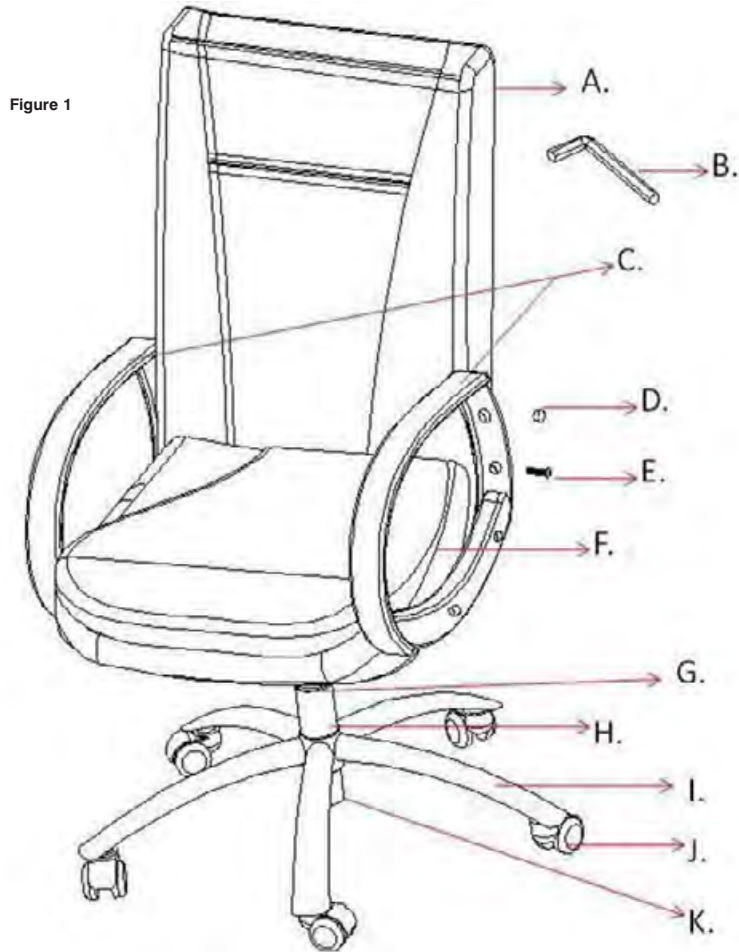
# Instructions

1. Place the star base face down and firmly press the rolling castors into place.

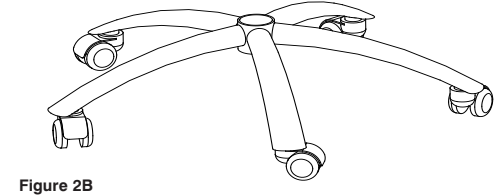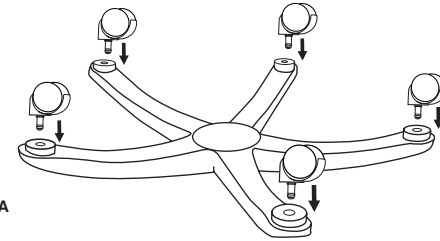2. Insert the pneumatic cylinder into the star base and cover with the plastic cone cover.

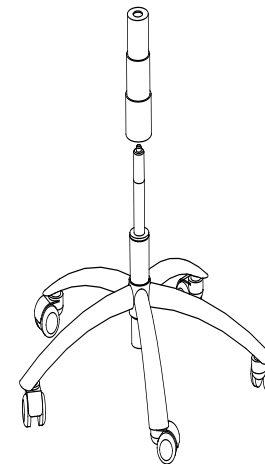3. Attach the seat base mechanism to the underside of the seat cushion using the included Allen wrench and 4 of the screws. Make sure the arrow on the seat base mechanism is pointing toward the front edge of the seat cushion (label on seat cushion).

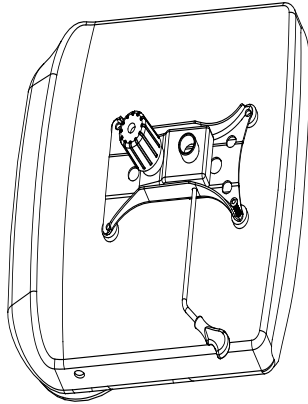
Figure 4

4. Place the seat cushion on top of the pneumatic cylinder. Make sure the cylinder is positioned properly by sitting on the seat cushion.

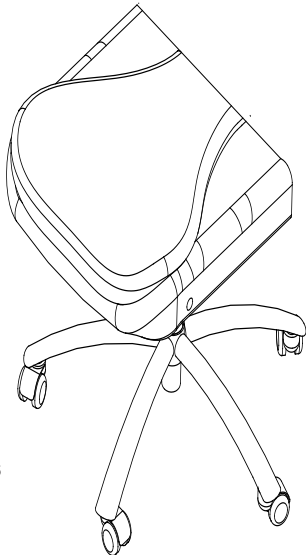
Figure 5

5. Attach the armrests (labeled as R and L for right and left sides) onto the seat cushion using the included Allen wrench and 4 of the screws (2 on each side). DO NOT TIGHTEN THE SCREWS COMPLETELY.

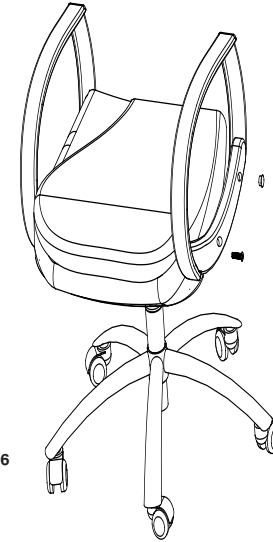
Figure 6

6. Connect remote arm wires to chair wires and feed into the opening in the side of the chair.

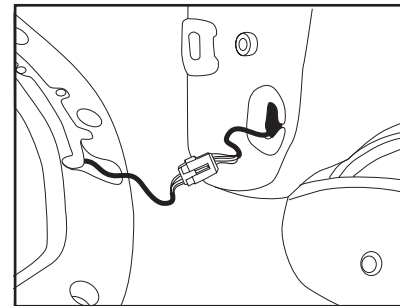
Figure 7A

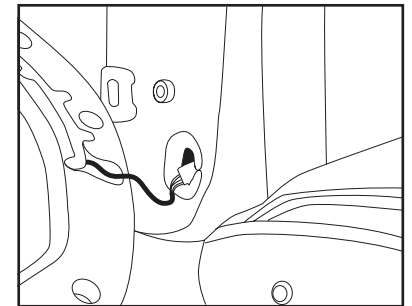
Figure 7B

7. Attach the seat back onto the armrests so it is properly aligned with the seat cushion using the included Allen wrench and the remaining 4 screws (2 on each side).

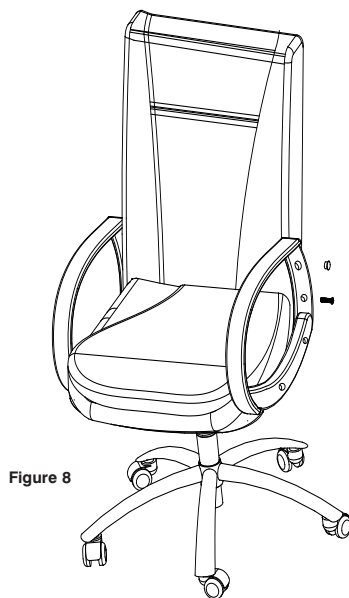
Figure 8

8. Once the armrests are properly aligned with both the seat cushion and seat back, completely tighten all 8 screws (4 on each side) and cover with provided screw covers.

**Note:**  
**Heat only works with massage function activated**

## Massage Cushion and Hand Control Instructions

### Power On/Off

To turn the massage function on, press the power button. The corresponding LED will illuminate to indicate power on. To turn the massage function off, press the power button again. As the chair is powering off, the power button LED will flash and the mechanism will automatically move to the bottom of the chair where it will “park” for the next use. No other buttons or functions will be available during the parking. Once the mechanism has parked, the power LED will turn off.

### Shiatsu Programs

Once the chair is powered on, you can select a massage program.

**Upper Back:** Press the button labeled “Upper Back” to select the massage program for that area. The corresponding LED will illuminate to indicate that program is engaged. To end the upper back program, press the button a second time. The LED will turn off.

**Lower Back:** Press the button labeled “Lower Back” to select the massage program for that area. The corresponding LED will illuminate to indicate that program is engaged. To end the lower back program, press the button a second time. The LED will turn off.

**Full Back:** Press the button labeled “Full Back” to select the massage program for that area. The corresponding LED will illuminate to indicate that program is engaged. To end the full back program, press the button a second time. The LED will turn off.

### Heat

To use heat during any massage program, press the “Heat” button. The corresponding LED will illuminate to indicate that the heat function is on. To end the heat function, press the heat button a second time. The LED will turn off. Please note that the heat function CANNOT be used without the massage function.

### Demo

To cycle through all the massage programs and heat, press the “Demo” button. The corresponding LED will be flashing to indicate the demo mode is functioning. Please note that if no other massage program is selected, the massage mechanism will park at the end of the demo cycle.

# Shiatsu Massage Chair

## Traveling Shiatsu Mechanism

### Massage Program Buttons

Direct massage to specific areas: lower back, upper back or full back. Press the button for your desired program and the corresponding L.E.D. will illuminate. To deselect an area(s), simply push the button again and the corresponding L.E.D. light will turn off.

### Power Button

To turn on the massage functions, first press the power button. The L.E.D. indicator will illuminate and confirm your selection. The massage functions can then be activated by pressing the appropriate buttons. To turn off the massage functions, simply press the button again and the L.E.D. indicator will turn off.

### Demo Button

For a brief demonstration of the features of the Traveling Shiatsu Office Chair, press the demo button. The unit will complete one cycle of each massage style. Once complete, the unit automatically turns itself off.

### Heat

Soothing heat is activated with this button. To activate the heat, first turn on the massage function. The L.E.D. light will glow to indicate the heat is on. To turn off heat, simply press the button again.

### Note:

As with all high-quality massage products, there is a 15-minute auto shut-off on the massaging mechanism.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Note: The manufacturer is not responsible for any radio or TV interference caused by unauthorized modifications to this equipment. Such modifications could void the user authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one of more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

## Product Care & Maintenance

### To Clean

- Disconnect the office chair from power source before cleaning.
- A soft, damp cloth may be used for the chair. When cleaning the controller use a dry cloth only.
- Dust may be removed from the backrest or armrest with a vacuum.
- Do not use abrasive cleaners, brushes, gasoline, kerosene, polish or paint thinner to clean.
- The cushions are not machine washable. Do not immerse in liquids to clean.

### To Store

- Avoid contact with sharp edges or pointed objects, which might damage/puncture the fabric.
- Do not store near heat or open flame. Do not leave this product exposed to direct sunlight for long periods of time. This may cause fading.

## Other General Care & Use Guidelines

- Foreign objects should not be inserted between the massage nodes or the unit housing.
- The mechanics of the product are specially designed and constructed to be maintenance free. No periodic lubrication or servicing is required.
- To retract cord into chair for storage or transport, unplug the chair from the wall. Loosely hold the cord and press the cord retraction button just above the cord on the chair. Gently guide the cord back into the chair while holding the retraction button.

### Tilt Lock:

To engage the tilt lock, make sure the chair is fully upright. Grasp the black handle located beneath the seat on the right side of the chair and push straight in (toward the chair). The tilt lock mechanism will click into place. The chair is now unable to recline. To disengage the tilt lock, make sure the chair is fully upright and pull the handle straight out (away from the chair). The chair is now able to recline.

### Height Adjust:

To raise the height of the chair, grasp the black handle located beneath the seat on the right side of the chair and lift up. Do not sit in the chair or put pressure on the seat while you are raising the handle. Release the handle when the chair is at the desired height. To lower the seat height, sit in the chair and grasp the black handle located beneath the seat on the right side of the chair and pull straight up. The chair will lower due to the weight on the chair. Release the handle when the chair reaches the desired height.

## Troubleshooting

<b>Problem</b>	<b>Probable Cause(s)</b>	<b>Suggested Solution</b>
<b>Massage action does not start</b>	No massage function was pressed.	Press one of the massage function buttons to start manual or programmed massage.
	The power cord is not plugged in.	Plug the power cord into a 120-volt AC outlet.
	The power button on the chair is not turned on.	Press the power button.
<b>Desired massage action does not start</b>	The correct massage function button was not pressed.	When pressing the massage function buttons, watch the LED's on the control. They will indicate the active function.
	The massage function button was pressed too quickly after the last button.	Wait several seconds and press the desired function again.
<b>Massage rollers stop while in action</b>	The self-contained safety device may stop the action if the chair has operated for an excessive period of time.	Unplug the chair and wait 30 minutes before restarting.
<b>Massage action does not stop after normal 15 minute auto shut off cycle</b>	Additional press of massage function button could have restarted another massage cycle.	Press the power button.

## Limited Warranty

In the unlikely event that your chair malfunctions or requires repair, please contact HoMedics Consumer Relations at 800-466-3342 or [cservice@homedics.com](mailto:cservice@homedics.com) between the hours of 8:30AM and 5:00PM EST for assistance. The customer should have the proof of purchase on hand. A trouble shooting process is completed by the customer and the consumer relations representative. Service will be scheduled or replacement parts sent if the chair is deemed to need warranty service. All warranty service requires approval and authorization by HoMedics. The determination of whether warranty service is necessary is made at HoMedics' sole discretion.

## LIMITED TWO YEAR WARRANTY

HoMedics sells its products with the intent that they are free of defects in manufacture and workmanship for a period of two years from the date of original purchase, except as noted below. HoMedics warrants that its products will be free of defects in material and workmanship under normal use and service. This warranty extends only to consumers and does not extend to Retailers.

To obtain warranty service on your HoMedics product, mail the product and your dated sales receipt (as proof of purchase), postpaid, to the following address:

HoMedics Consumer Relations  
Service Center Dept. 168  
3000 Pontiac Trail  
Commerce Township, MI 48390

No COD's will be accepted.

HoMedics does not authorize anyone, including, but not limited to, Retailers, the subsequent consumer purchaser of the product from a Retailer or remote purchasers, to obligate HoMedics in any way beyond the terms set forth herein. This warranty does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorized accessory; alteration to the product; improper installation; unauthorized repairs or modifications; improper use of electrical/power supply; loss of power; dropped product; malfunction or damage of an operating part from failure to provide manufacturer's recommended maintenance; transportation damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of HoMedics.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adoption to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. THERE SHALL BE NO OTHER WARRANTIES EXPRESS OR IMPLIED INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS OR ANY OTHER OBLIGATION ON THE PART OF THE COMPANY WITH RESPECT TO PRODUCTS COVERED BY THIS WARRANTY. HOMEDICS SHALL HAVE NO LIABILITY FOR ANY INCIDENTAL, CONSEQUENTIAL OR SPECIAL DAMAGES. IN NO EVENT SHALL THIS WARRANTY REQUIRE MORE THAN THE REPAIR OR REPLACEMENT OF ANY PART OR PARTS WHICH ARE FOUND TO BE DEFECTIVE WITHIN THE EFFECTIVE PERIOD OF THE WARRANTY. NO REFUNDS WILL BE GIVEN. IF REPLACEMENT PARTS FOR DEFECTIVE MATERIALS ARE NOT AVAILABLE, HOMEDICS RESERVES THE RIGHT TO MAKE PRODUCT SUBSTITUTIONS IN LIEU OF REPAIR OR REPLACEMENT.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or sales of such products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate as to any products or parts thereof which are repaired, replaced, altered, or modified, without the prior express and written consent of HoMedics.

This warranty provides you with specific legal rights. You may have additional rights which may vary from country to country. Because of individual country regulations, some of the above limitations and exclusions may not apply to you.

For more information regarding our product line in the USA, please visit: [www.homedics.com](http://www.homedics.com)

**HoMedics Service Center**  
**1.800.466.3342**  
**8:30 – 5:00 p.m. (EST) M – F**

**Mail To:**  
**HoMedics**  
**Service Center Dept. 168**  
**3000 Pontiac Trail**  
**Commerce Township, MI**  
**48390**

**Email:**  
**[cservice@homedics.com](mailto:cservice@homedics.com)**

©2010 HoMedics, Inc. and its affiliated companies, all rights reserved. HoMedics® is a registered trademark of HoMedics Inc. and its affiliated companies. All rights reserved.  
1B-OCTS200

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