
Ozeri PRECISION

Digital Bathroom Scale

User Manual

Model: ZB18

(Model # may vary by color)

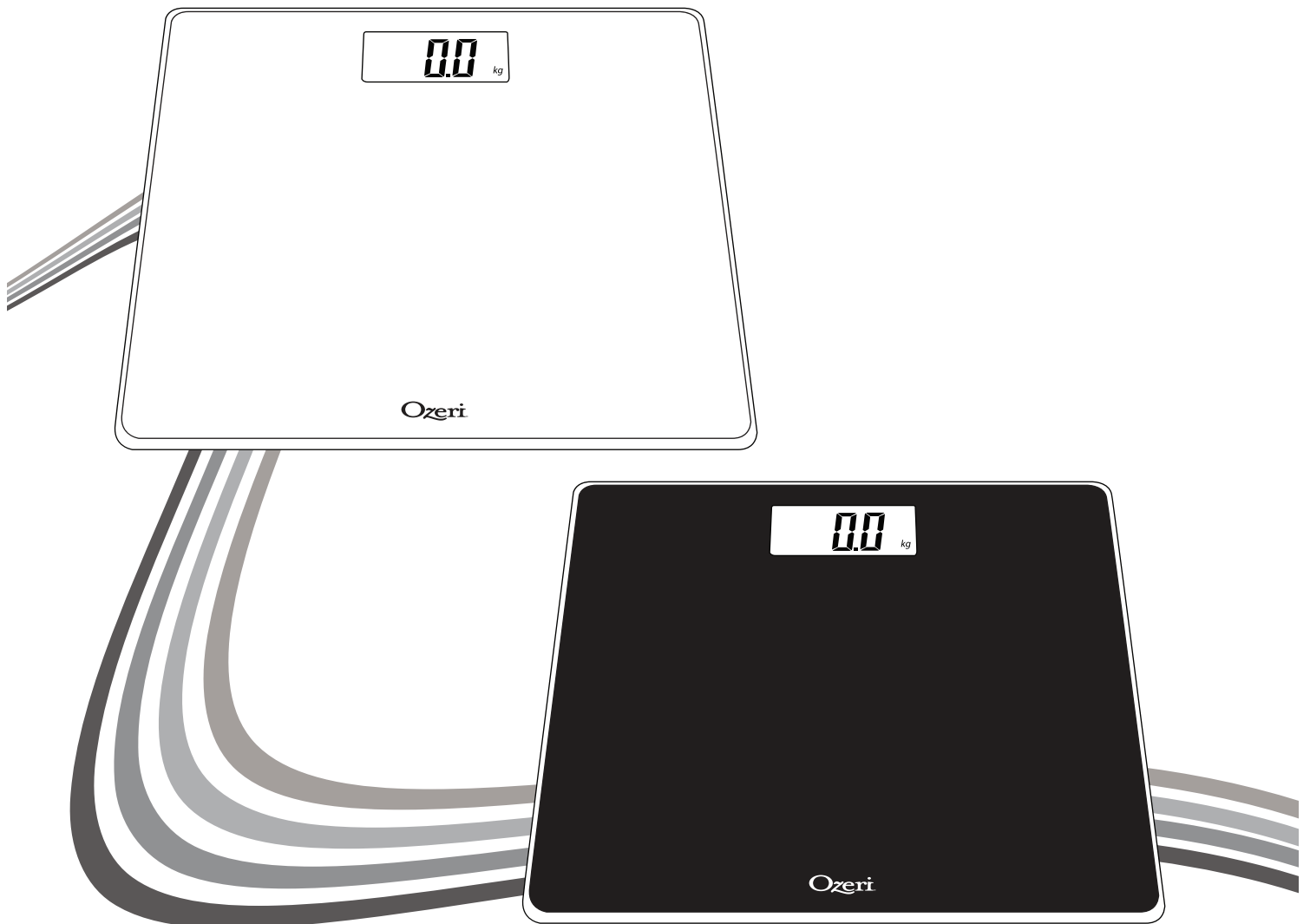Congratulations

Congratulations on your purchase of the Ozeri Precision Digital Bathroom Scale. Please read through this manual before first using your new scale, and keep this manual handy for future reference.

Warnings


- Your Ozeri Precision Digital Bathroom Scale is a finely tuned electronic instrument. One of the hardest things to weigh is the human body because there is always subtle motion. To ensure reliable readings, always try to place your feet on the scale in the same place each time and remain still.
- Your scale will auto-calibrate once before first use and after each time the scale detects that it has been picked up or moved to another location. As part of this calibration process, please disregard the first weight reading. Step off the scale and allow it to turn off. This will complete the calibration process.
- Note that your body weight will vary based on hydration, clothing and other factors. This scale rounds to the nearest 100 grams. When slight variances occur, it is recommended to average your weight over time.
- Treat your scale with care. Do not drop it or jump on it. While the scale is made out of durable materials, misuse or abuse may render the electronic sensors inoperable.
- The maximum weight capacity of this scale is 440 lbs. To avoid permanent damage, do not attempt to exceed the weight limit and do not store any items on this scale when not in use.
- Do not disassemble this product. Other than replacing the batteries, it contains no user serviceable parts.

- The product is intended for home / consumer use only.
- Clean after use with a slightly dampened cloth. Do not use solvents or immerse the product in water. Thoroughly dry the scale after cleaning. Moisture can damage the scale's sensors and electronics.
- Do not use the scale on a carpet or uneven surface.
- Do not store the scale where you store cleaning chemicals. The vapors from some household products may affect the electronic components of your scale. Do not store the scale on its side.
- Store your scale in a clean, dry location at room temperature. Dust, dirt, and moisture from humidity can accumulate on weighing sensors and electronics, causing inaccuracy or malfunction.
- To prevent battery drainage, do not store anything on the scale.
- Do not dispose of the batteries in fire. Batteries may explode or leak.

Specifications

- Tempered Glass 11.8" x 11.8" Platform
- Weight Units: lbs / kg / st
- Maximum Capacity: 440 lbs / 200 kg
- Minimum Weight Capacity: 13.2 lbs / 6 kg
- Weight Graduation: 0.2 lbs / 0.1 kg
- Power Supply: 2 AAA batteries (included)
- Low Battery Indicator: "Lo"
- Auto On / Off: Powers on when you step on the scale; powers off after 10 seconds of inactivity
- Equipped with 4 high precision sensors
- Tare button for weighing of infants, pets, or luggage

Battery

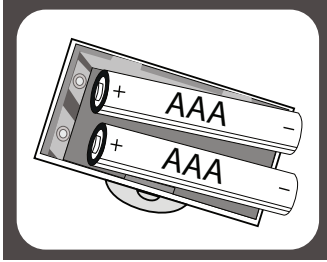Insert the included 2 AAA batteries into the battery compartment located underneath the scale. Ensure +/- terminals are correctly oriented. When replacing batteries, please use a small tool to help remove the batteries from the battery compartment if necessary.

Using the Scale

Note: Your scale will auto-calibrate once before first use and after each time the scale detects that it has been picked up or moved to another location. As part of this calibration process, please disregard the first weight, step off the scale and allow it to turn off. This will conclude the calibration process.

- Once calibrated, place your scale on an even hard flat surface. Do not use your scale on carpets or soft uneven surfaces. Ensure your scale is dry prior to use. **WARNING: NEVER** step onto a wet scale and always ensure solid footing to avoid injury.

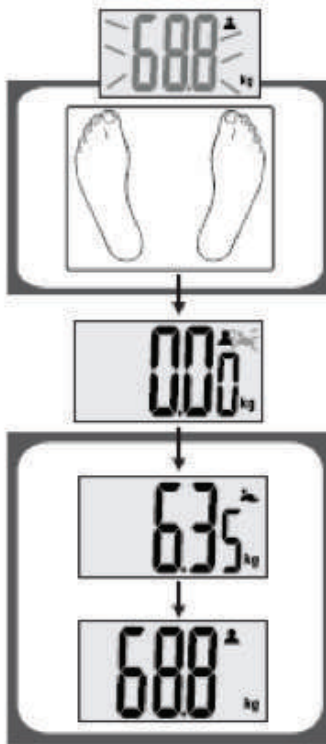
- Step onto the scale. Your scale features InstaOn™ technology that requires no tapping to activate the scale. One or more weight readings will display on the screen while your scale calculates your weight. When complete, your final weight calculation will remain on the screen.
- The scale will automatically turn off after approx 10 seconds of inactivity.

Tare Function

This scale features a Tare function that allows you to weigh your infant, pet or luggage piece by carrying it while standing on the scale.

To activate / deactivate the Tare function, first ensure the scale is on by briefly stepping onto it. Then press and hold the Unit button located underneath the scale for 3 seconds. When the Tare function is activated, the “” icon will appear on the top right side of the screen.

Step onto scale without carrying anything. One or more weight readings will display on the screen while your scale calculates your weight. When complete, your final weight calculation will flash three times and return to zero. The infant “” icon will then appear and start flashing on the screen.

The scale will wait 30 seconds for you to safely pick up and hold your infant (or other object) and step onto the scale. One or more weight readings will display on the screen while the scale calculates the infant’s weight. When complete, the infant’s (or other object’s) final weight will flash three times and remain on the screen for 10 seconds. The screen will then alternate back to show your weight again after showing the infant’s weight, and then power off.

Unit of Measurement

Your scale is set to calculate weights in pounds (lbs). To change the unit of measurement between pounds, kilograms and British stones, simply use the Unit switch located underneath the scale.

Troubleshooting Guide

If you experience any difficulties with your scale, please ensure the following:

- Double check that the batteries are correctly inserted according to their +/- polarity.
- Ensure the scale is on a flat, level floor and not touching a wall or other object.

- If the Low Battery indicator appears on the screen (Lo), replace battery with 2 new AAA batteries of the same type.
- If the “Err” appears on the screen, the maximum weight of 440 lbs has been exceeded.

WARRANTY INFORMATION

Your “Ozeri Scale” is warranted to the original owner for 1 year from the date of purchase against defects in material and workmanship under normal use and service. Should your scale prove defective within 1 year from the date of purchase, return the defective part or unit, freight prepaid along with an explanation of the claim to “Ozeri Scale.” Please pack your scale carefully to avoid damages in transit. Under this warranty, Ozeri will repair or replace any parts found to be defective due to a manufacturer’s defect. This warranty extends only to personal use and does not extend to any product that has been used for commercial use, rental use, or use which is not intended. There are no warranties other than that expressly set forth herein. This warranty is not transferable. Ozeri is not responsible in any way for damages, losses or inconveniences caused by equipment failure or by user negligence, abuse, or use not in accordance with the “User Guide” or any additional safety, use, or warnings included in the product materials. Further, Ozeri is not liable for incidental or consequential damages of any nature resulting from the use of this product and any liability shall not exceed the purchase price of the product. This Warranty gives you specific legal rights. You may also have other rights which vary from state to state.

THIS WARRANTY DOES NOT APPLY TO THE FOLLOWING

- Damage due to accident, misuse, abuse, alterations or vandalism.
- Improper or inadequate maintenance.
- Unauthorized modification or commercial use.
- Damage in return transit.
- Unsupervised use by children under 18 years of age.


Ozeri Customer Service:
Phone: 1-877-299-1296
Email: support@ozeri.com

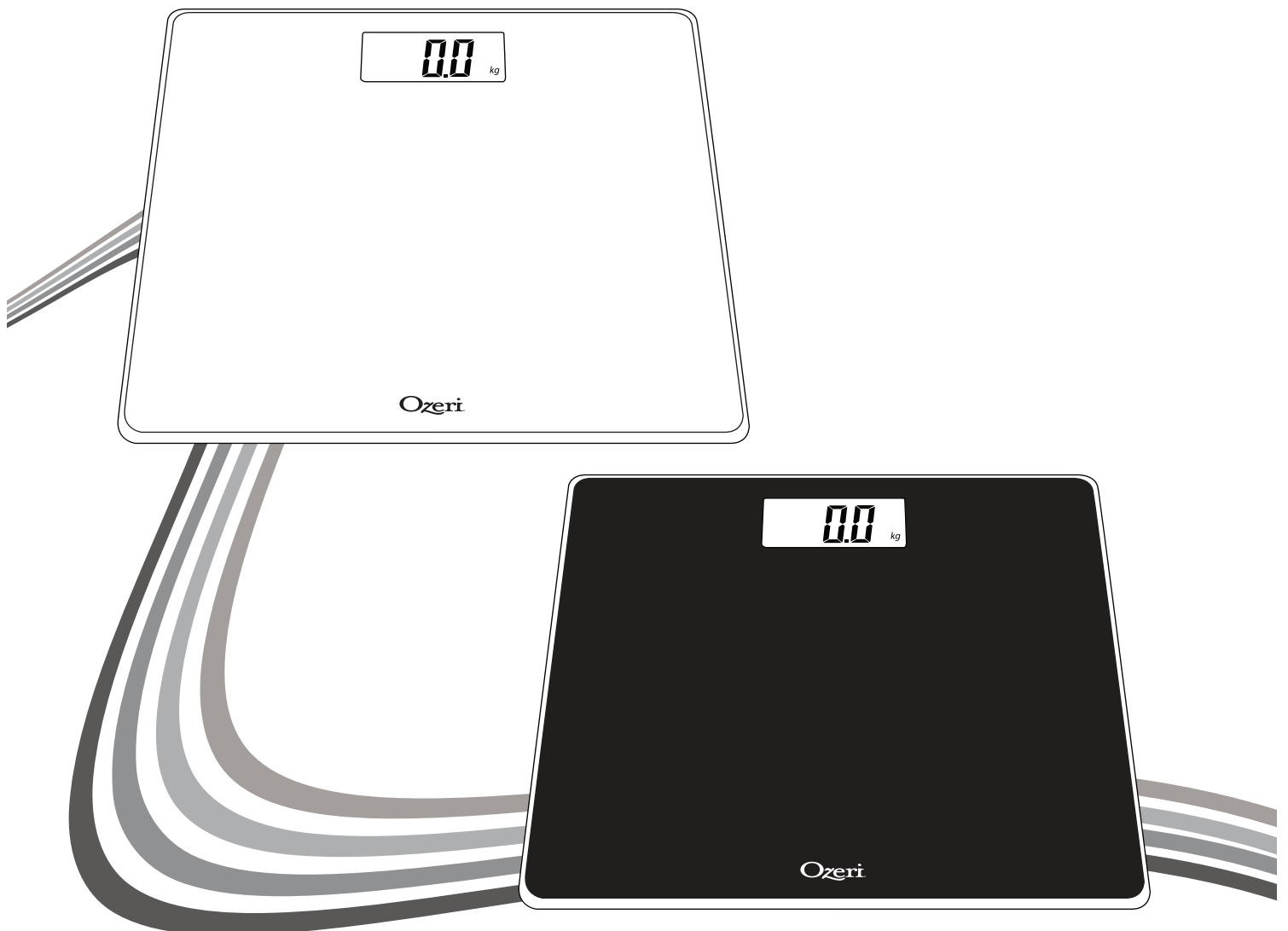(Ozeri[®] and its subsidiaries assume no liability for damage caused by the use of the Ozeri Scale other than for its intended use or as instructed above and in the User Manual.)

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